

pOET6 BacMam Transfer Plasmid

Product Data Sheet

Product information for MIR 6152

SDS and Certificate of Analysis available at mirusbio.com/6152



SPECIFICATIONS

Storage	Store pOET6 BacMam Transfer Plasmid DNA at -20°C .
Product Guarantee	1 year from the date of purchase, when properly stored and handled.
Concentration	20 μl plasmid DNA at 500 ng/ μl , sterile filtered in TE Buffer

Plasmid Description: pOET6 BacMam is a baculovirus transfer vector designed to facilitate the expression of foreign genes in mammalian cells under the cytomegalovirus (CMV) immediate early gene promoter. Vector features include an ampicillin resistance gene for selection in *E. coli*, a Col E1 origin of replication, a CMV promoter to drive constitutive foreign gene expression, and a multiple cloning site (MCS) containing 11 unique restriction sites for inserting the gene in the correct orientation. Following transfection into insect cells, the AcMNPV and ORF1629 sequences flanking the expression cassette allow recombination into the baculovirus genome which then replicates to produce baculovirus that can be harvested directly from the cell culture medium. Purified baculovirus can be used to transduce a variety of mammalian cells in a titer-dependent manner. pOET6 BacMam is compatible with the *flashBAC*[™] System or any baculovirus expression system that utilizes homologous recombination at the polyhedrin locus in insect cells.

Plasmid Features:

AcMNPV: 329 - 1156

CMV (promoter): 1380 - 2188


Multiple Cloning Site (MCS): 2252 - 2303

SV40 poly(A) Signal: 2313 - 2434

ORF1629: 2580 - 3921

Col E1 Replication Origin: 4169 - 4788

Amp Resistance Marker: 4943 - 5803

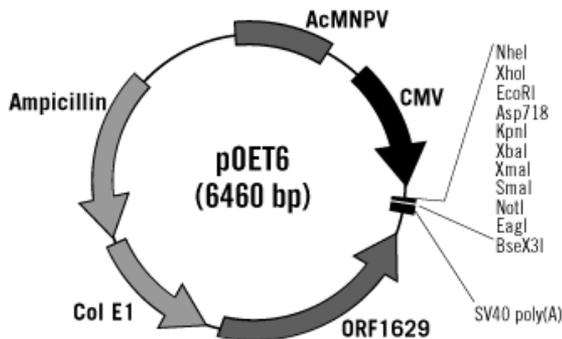

Additional Information: All pOET vectors can be propagated in DH5 α or other general purpose *E. coli* strains. For more information on baculovirus expression using the *flashBAC*[™] System, pOET transfer plasmids, or optimized insect cell transfections using *TransIT*[®]-Insect Transfection Reagent, please visit www.mirusbio.com.

For Research Use Only

Mirus Bio LLC

www.mirusbio.com | techsupport@mirusbio.com | Toll Free (U.S.): 844.MIRUSBIO | Direct: +1.608.441.2852

► NOTES



SDS and Certificate of Analysis available at [mirusbio.com/6152](https://www.mirusbio.com/6152)

flashBAC™ Systems and pOET vectors are sold by Mirus Bio through partnership with Oxford Expression Technologies, Oxford, UK.

©1996-2018 All rights reserved. Mirus Bio LLC. All trademarks are the property of their respective owners.
For terms and conditions, visit www.mirusbio.com

Rev.A 1216

Mirus Bio LLC

www.mirusbio.com | techsupport@mirusbio.com | Toll Free (U.S.): 844.MIRUSBIO | Direct: +1.608.441.2852